


CO-EFFICIENT OF DISCHARGE APPARATUS MAKE – LABTEK, ORIGIN – INDIA

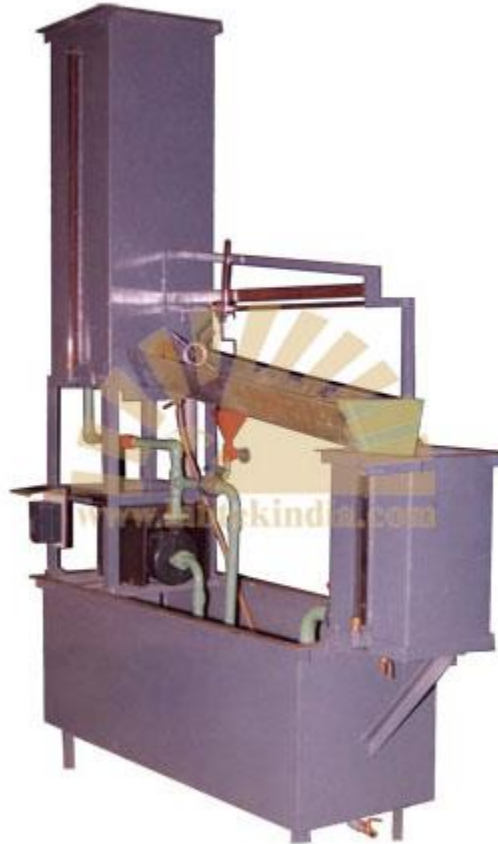

Supply Tank:

M.S Powder Coating/M.S.FRP Lining/S.S Unit of suitable size with a fixture arrangement to mount 2 replaceable orifices. A scale and Sliding arrangement for measurement of X and Y ordinate to find the velocity traverse of the water jet from orifice and provided with a micrometer contraction gauge for the measurement of contraction of water jet at the throat of the orifice.

Orifices :

Two Orifices made of gunmetal