


DEVAL ABRASION TESTING MACHINE MODEL NO SL -BT -011, MAKE - LABTEK, ORIGIN - INDIA

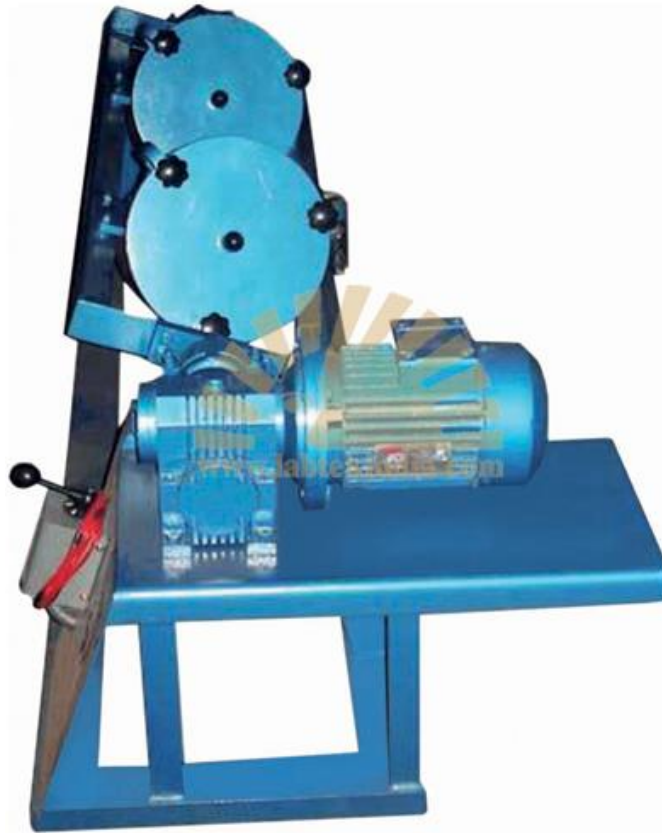

ASTM D2-33 & IS 2386 (PART IV)

For The Determination of Resistance of Aggregates to Wear by Abrasion.

Specification :

It consists of two hollow cylinders closed at one end and provided with Fitting covers at the either end. These cylinders are mounted on a shaft at an angle of 30° with the axis of rotation of the shaft. The shaft rotates at 30-33 RPM. Through a reduction gear operated by a motor and is provided with a revolution counter. Complete with Abrasive Charge consisting of 12 Nos. Hardened steel Balls of 48mm dia. Suitable for operation on 440 Volts, Three Phase, 50 Cycles, A.C. supply.