


DUCTILITY TESTING APPARATUS

MODEL NO SL -BT -015, MAKE - LABTEK, ORIGIN - INDIA

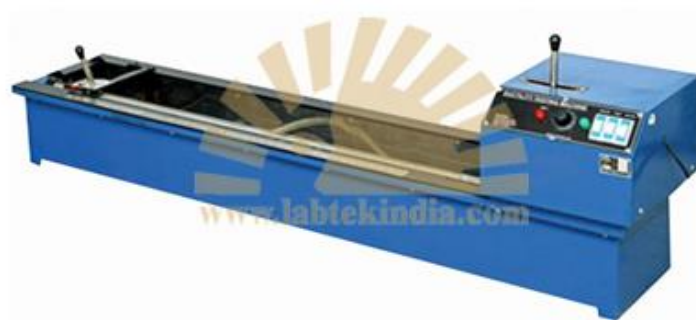

IS 1208-1058, ASTM D 113, IP32, 55, AASHTO T 51

Designed to test three specimens simultaneously. The machine consists of a carriage moving over a lead screw. An electric motor driven reduction gear unit ensures smooth constant speed and continuous operation. The entire assembly is mounted with a stainless steel lined water bath completely encased in metal bound hardwood. It is equipped with an electric pump circulator and heater. The temperature is controlled thermostatically. Two rates of travel i.e. 5 cm/min and 1 cm/min are provided. Suitable for operation on 230 V, 50 Hz, Single Phase, A.C. supply.